




RELIEVE YOUR PAIN WITH MAGNET THERAPY Relieve your pain with our **Bio-magnets 1500G**

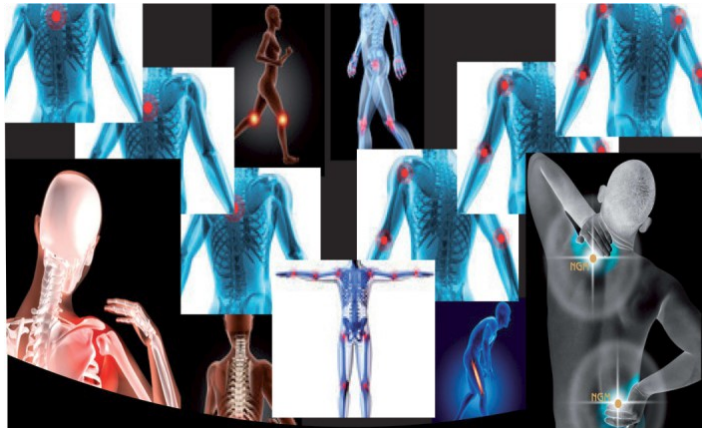

Shoulder condition

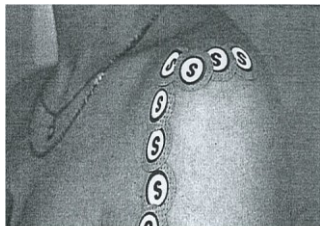

Unstable ankle

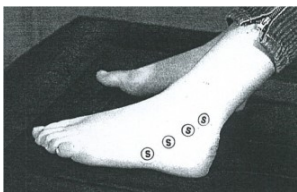

Fracture of the midfoot



Prosthesis on both knees. One year later:
algodystrophy, pain, ...
Impossibility to climb stairs.
Placement of magnets on each knee.
Ability to climb stairs has returned.

Broken arm:


2 days after surgery



Day +2: placement of 12 magnets
against the pain and to speed up
calcification



Arthrosis of the knee

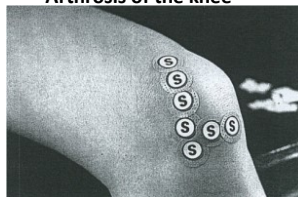

Major applications and action:

The magnetic fields help suppress joint, muscle and rheumatic pain. They have a calming and soothing action on a great diversity of pains, but also on control of stress . Placed on the skin with the:

- Northside: relaxing and pain relieving effect
- Southside: anti-inflammatory action.

Proven efficiency – clinically tested:

The totality of testing in many countries with the bio-magnets 1500G show an efficiency of 85% for pain relief.
Completely harmless: no side effects.

1– How to apply:

Determine the exact location of the pain points through palpation at rest and during movement or with the foldable Lecher antenna cereb ©. Then put the bio-magnet 1500 g on an adhesive and put on the skin. Verify whether there are other pain points and continue as before.

2– Advantages of the bio-magnets 1500G:


- These miniaturised bio-magnets 1500 g with permanent magnetization and constant intensity are reusable for many years.
- the North and South Pole are referred to
- easy to use on any part of the body

3– Precautions: DO NOT use:

- in open wounds, in persons allergic to metal
- in persons with a pacemaker
- in pregnant women and young children
- for external use only

Specifications bio-magnets 1500G:

Southside on the skin : anti-inflammatory
(=side with blue indication)

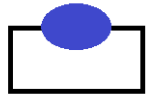

Northside on the skin: relaxing

1– Characteristics:

The strong bio-magnets are offered in the form of magnetized pastilles of reduced size (5 mm/1 mm) in order to facilitate the application and this regardless of the zone to be treated on the body.

- consisting of rare-earth materials they are considered permanent magnetic as they don't lose more than 1/1000 of a % of their magnetic flux per year:
- the magnet is plated with half a micron of gold to avoid oxidation and any other cutaneous intolerance, as well as the risk of allergy to metal alloys.

Place, in case of itchiness an additional adhesive tape between the skin and the magnet.

The South side of the magnet is recognizable by the blue indication on one of the sides. The magnet is put on an adhesive tape and then glued at the height of the pain point.

Order:

5 magnets for 29 Euros

10 magnets for 49 Euros

Your guarantee: satisfied or money back

WWW.CEREB.BE

Mobile number: 00 32 497 467 699

e-mail: CEREB@skynet.be for enquiries

Ask for the documentation :

Bio-magnets 1500G and acupuncture

MAGNETS FOR YOUR HEALTH

After each number of the point the Chinese name is mentioned and the exact reference number from traditional acupuncture (LU9 = 9th point of the Lung Meridian). The magnetized pastille (several square millimeters) does not need to be applied to the acupuncture point as with the same precision as a needle. The following schemes are therefore sufficient. For more precise localization reference is made to The Concise Book of Acupoints by John R. Cross.

0 Zhangmen (LR13): STOMACH PROBLEMS

1 Taiyang (OM 2): MIGRAINE

2 Shufu (KI 27): BRONCHIAL ASTHMA

3 Yuzhong (KI 26) + Shencang (KI 25): BRONCHIAL ASTHMA ,
COLD, SOAR THROAT, COUGH

4 Jiuwei (CV 15): HEADACHE, MIGRAINE, NERVOUS DISORDER

5 Shousanli (LI 10): SHOULDER PAINS

6 Quchi (LI 11): LARYNGITIS, LUMBAGO, MEMORY LOSS

7 Shenjue (CV 8): INTESTINALE PARESIS, STOMACH COMPLAINTS, NAUSEA, CONSTIPATION, BLOATING

8 Qihai (CV 6): INSOMNIA

9 Zhongji (CV 3): COLIBACILLOSIS

10 Dubi (ST 35): HIP COMPLAINTS, HAEMORRHOIDS

11 Jianzhongshu (SI 15): SHOULDERPAIN, HEADACHE

12 Jianwaishu (SI 14): IDEM TO ABOVE

13 Shenshu (BL 23): LUMBAGO, REUMATIC COMPLAINTS

14 Chanqiang (GV 1): HAEMORRHOIDS

15 Ququan (LR 8): ASTHMA, FRIGIDITY, HEPATIC INSUFFICIENCY, MIGRAINE

16 Weizhong (BL 40): ANEMIA

17 Tianjing (TE 10): AGITATION, NERVOUSNESS

18 Wenliu (LI 7): IDEM AS ABOVE

19 Waiguan (TE 5): ELBOW AND SHOULDER COMPLAINTS

20 Taiyuan (LU 9): ABCESS, TONSILLITIS, BRONCHITIS

21 Yangchi (TE 4): ANEMIA, DIABETES

22 Yangxi (LI 5): LUMBAGO

23 Zhongzhu (TE 3): HEADACHE

24 Hegu (LI 4): BRONCHITIS, CONSTIPATION, MIGRAINE

25 Sanjian (LI 3): ACNE

26 Erjian (LI 2): ACNE, TONSILLITIS, LARYNGITIS

27 Shaohai (HT 3): FEAR

28 Chize (LU 5): ASTHMA

29 Shaoze (SI 1): CONSTIPATION

30 Guanchong (TE 1): MIGRAINE

31 Shangyang (LI 1): BRONCHITIS, ACNE, TOOTHACHE

32 Shenmen (HT 7): ANXIETY

33 Jingqu (LU 8): ASTHMA, BRONCHITIS

34 Shaofu (HT 8): ASTHMA, HYPERTENSION, MEMORY LOSS

35 Shaoshang (LU 11): TONSILLITIS, LARYNGITIS

36 Shaochong (HT 9): ANXIETY, HYPERTENSION, STAGE FRIGHT

37 Taixi (KI 3): ABCESS

38 Zhaohai (KI 6): ANKLE PAIN

39 Rangu (KI 2): DIABETES

40 Dadu (SP 2): ANOREXIA, DIABETES

41 Yinbai (SP 1): ANEMIA

42 Yangfu (GB 38): HAEMORRHOIDS

43 Shenmai (BL 62): ANKLE PAIN

44 Zhongfeng (LR 4): ANOREXIA, STOMACH COMPLAINTS

45 Qiuxu (GB 40): HIP COMPLAINTS

46 Shangqiu (SP 5): HAEMORRHOIDS

47 Xiangqu (ST 43): HYPERTENSION

48 Taichong (LR 3): FRIGIDITY, IMPOTENCE, HEPATIC INSUFFICIENCY

49 Qiaoyin van de voet (GB 44): ABCESS

50 Tonggu (BL 66): HIP COMPLAINTS

51 Shugu (BL 65): ABCESS, HIP COMPLAINTS

52 Jinggu (BL 64): ANOREXIA

53 Yongquan (KI 1): MEMORY

54 Special point for ASTHENIA AND GENERAL WEAKNESS

55 Tongli (HT 5): STAGE FRIGHT

